



### AccuDiag™ NSE ELISA KIT

REF 6334-16

IVD See External Label 8°C 96 Tests

Enzyme Immunoassay for the Quantitative Measurement of NEURON-SPECIFIC ENOLASE (NSE) in Human Serum.

NSE ELISA	
Method	Enzyme Linked Immunosorbent Assay
Principle	Sandwich Complex
Detection Range	0-200 ng/ml
Sample	25 ul
Incubation Time	80 minutes
Shelf Life	12-14 Months from the manufacturing date

### PRODUCT FEATURES

- Very easy to use with little training
- Highly specific and consistent Assay
- Provides accurate results quickly
- Reading of results both visually and as absorbance data

### INTENDED USE

The Neuron-Specific Enolase (NSE) enzyme linked immunosorbent assay (ELISA) provides quantitative measurement of human NSE in serum to aid in the clinical evaluation of patients suspected of having small cell lung cancer, and other related diseases.

### SIGNIFICANCE AND SUMMARY

The glycolytic enzyme enolase (2-phosph-D-glycerate hydrolyase) exists as several dimeric isoenzymes ( $\alpha\alpha$ ,  $\alpha\beta$ ,  $\alpha\gamma$  and  $\gamma\gamma$ ) composed of three distinct subunits,  $\alpha$ ,  $\beta$ , and  $\gamma$ . Three isoenzymes are found in human brain:  $\alpha\alpha$ ,  $\alpha\gamma$  and  $\gamma\gamma$ . The  $\alpha\gamma$  and  $\gamma\gamma$ -enolase isoenzymes are also known as neuron-specific enolase (NSE) as these isoenzymes initially were detected in neurons and neuronendocrine cells. The NSE levels are low in health and benign subjects.

Elevated levels are commonly found in patients with malignant tumors with neuronendocrine differentiation, especially small cell lung cancer and neuroblastoma.

Lung cancer is one of the most spread cancer forms with incidences about 50-100 per 100,000 population. Approximately 20% of the lung cancer is small cell lung cancer. Patients with small cell lung cancer show various proportions of  $\alpha\gamma$  and  $\gamma\gamma$  isoenzyme. The determination of NSE should detect and  $\gamma$  isoforms with the same sensitivity (1). The antibodies for this particular assay are specific for the  $\gamma$ -subunit without cross reactivity with  $\alpha$  or  $\beta$  subunits(1). NSE are reported to be useful diagnostic marker for lung cancer(2), neuroblastoma(3), melanoma(4), seminoma(5) and in injury of central nervous system(6). In addition to the above, NSE can be a valuable tool in following-up the effect of chemotherapy of small cell lung cancer, in prognostic evaluation of patients with small cell lung cancer, and in differential diagnosis between cell lung cancer and non-small cell lung cancer.

### ASSAY PRINCIPLE

The NSE Quantitative Test Kit is based on a solid phase enzyme-linked immunosorbent assay. The assay system utilizes one monoclonal anti- $\gamma$ NSE antibody for solid phase (microtiter wells) immobilization and another monoclonal anti- $\gamma$ NSE antibody in the antibody-enzyme (horseradish peroxidase) conjugate solution. The standards and test specimen (serum) are added to the antibody coated microtiter wells. During the incubation, specific NSE bound to anti-NSE antibody on the wells. Unbound NSE antigen is removed by washing the wells with wash buffer. Enzyme conjugate is then added to each well. After another incubation, unbound enzyme conjugate is washed off and the amount of bound peroxidase is proportional to the concentration of the NSE present in each sample. Upon addition of the substrate and chromogen, the intensity of blue color will develop in proportion to the concentration of NSE antigen in the samples.

### SPECIMEN COLLECTION AND PREPARATION

1. Blood should be drawn using standard venipuncture techniques and the serum should be separated from the red blood cells as soon as practical.
2. Plasma samples collected in tubes containing EDTA, heparin, or oxalate may interfere with test procedures and should be avoided.
3. For prolong storage of whole blood can cause release of NSE from the blood cells. If sera can not be assayed immediately, they can be stored at 2-8°C for one day or frozen at -20°C for up to 30 days prior to assay. Sample should not be refrozen. Repeated freezing and thawing is not recommended. Do not store in self-defrosting freezer.
4. Do not use hyperlipemic, hemolyzed, Plasma is not recommended since significant amounts of NSE can be released from the thrombocytes.
5. Avoid turbid and contaminated samples.

### MATERIALS AND COMPONENTS

#### Materials provided with the test kit

1. Monoclonal anti-NSE antibody coated microtiter plate with 96 wells.
2. Sample diluent (12 ml)
3. Enzyme conjugate reagent, 12 ml.
4. Lyophilized NSE reference standards containing; 0, 5, 15, 40, 100, and 200 ng/ml of NSE. 1 set.
5. 50 X Wash Buffer Concentrate, 15 ml.
6. TMB Substrate, 12 ml.
7. Stop Solution 12ml.
8. Control Set (optional)

Diagnostic Automation/Cortez Diagnostics, Inc.

21250 Califa St, Suite 102 and 116, Woodland Hills, CA 91367 USA Phone: 818-591-3030, Fax: 818-591-8383

Email: [onestep@rapidtest.com](mailto:onestep@rapidtest.com) Website: [www.rapidtest.com](http://www.rapidtest.com)



### Materials required but not provided

1. Precision pipettes and tips, 0.025ml, 0.05 ml, 0.10 ml, and 1.0 ml.
2. Disposable pipette tips.
3. Distilled water.
4. Glass tubes or flasks to prepare wash buffer.
5. Vortex mixer.
6. Absorbent paper or paper towel.
7. Microtiter plate reader.
8. Graph paper.

If a serum specimen contains greater than 200 ng/ml of NSE the sample must be diluted with sample diluent and re-assayed as described in the assay procedure.

### RESULTS

Calculate the mean absorbance value for each set of NSE reference standards, specimens and controls. Construct a standard curve by plotting the mean absorbance obtained from each reference standard against its concentration in ng per ml on linear graph paper, with absorbance values on the vertical or Y axis and concentrations on the horizontal or X axis. Use the mean absorbance values for each specimen to determine the corresponding concentration of NSE in ng per ml from the standard curve. Any diluted specimens must be corrected by the appropriate dilution factor.

### REAGENT PREPARATION


1. All reagents should be brought to room temperature (18-22°C) and mixed by gently inverting or swirling prior to use. Do NOT induce foaming.
2. Add 0.5 ml of distilled water to reconstitute the lyophilized standards. Allow the reconstituted materials to stand for at least 20 minutes. Mix gently. The reconstituted standards should be stored sealed at 2-8°C
3. Dilute 1 volume of Wash Buffer Concentrate (50x) with 49 volumes of distilled water. For example, Dilute 15 ml of Wash Buffer Concentrate (50x) into 735 ml of distilled water to prepare 750 ml of washing buffer (1x). Mix well before use.

### EXAMPLE OF STANDARD CURVE

Results of a typical standard run with optical density reading at 450nm shown in the Y axis against NSE concentrations shown in the X axis

NSE Conc (ng/ml)	Absorbance (450nm)
0	0.009
5	0.176
15	0.376
40	0.835
100	1.782
200	2.978

This standard curve is for the purpose of illustration only, and should not be used to calculate unknowns. Each user should obtain his or her own standard curve and data.

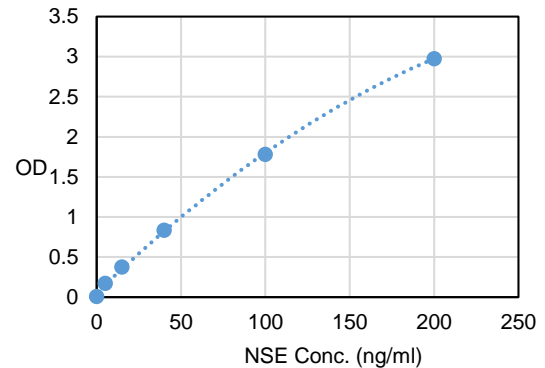

### EXPECTED VALUES AND SENSITIVITY

### ASSAY PROCEDURE

1. Secure the desired number of coated wells in the holder.
2. Dispense 25µl of standard, specimens, and controls into appropriate wells.
3. Dispense 100µl of Sample diluent into each well.
4. Thoroughly mix for 10 seconds. It is very important to have complete mixing in this step.
5. Incubate at room temperature (18-22°C) for 30 minutes.
6. Remove the incubation mixture by flicking plate content into a waste container.
7. Rinse and flick the microtiter wells 5 times with washing buffer (1X).
8. Strike the wells sharply onto absorbent paper or paper towels to remove all residual water droplets.
9. Dispense 100µl of Enzyme Conjugate Reagent into each well. Gently mix for 5 seconds.
10. Incubate at room temperature for 30 minutes.
11. Remove the incubation mixture by flicking plate contents into a waste container.
12. Rinse and flick the microtiter wells 5 times with washing buffer (1X).
13. Strike the wells sharply onto absorbent paper to remove residual water droplets.
14. Dispense 100µl TMB substrate into each well. Gentle mix for 5 seconds.
15. Incubate at room temperature for 20 minutes.
16. Stop the reaction by adding 100µl of stop solution into each well.
17. Gently mix for 30 seconds to make sure that the blue color changes to yellow color completely.
18. Read optical density at 450nm with a microtiter reader within 15 minutes.

#### Important Note:

1. The wash procedure is critical. Insufficient washing will result in poor precision and falsely elevated absorbance readings.
2. It is recommended that no more than 32 wells be used for each assay run if manual pipetting is used since pipetting of all standards, specimens and controls should be completed within 5 minutes. A full plate of 96 wells may be used if automated pipetting is available.
3. Duplication of all standards and specimens, although not required, is recommended.

1. It is recommended that each laboratory should determine its own normal and abnormal ranges as to account for its environmental factors such as diet, climate etc.
2. A clinical study of the NSE Quantitative kit was conducted and results are summarized as follows: Nearly all the individuals have NSE values below 15 ng/ml (95<sup>th</sup> percentile).
3. The expect ranges are representative only, and do not necessarily reflect the ranges that will be observed in a particular clinical laboratory.



### LIMITATIONS

There are some limitation of the assay:

1. As with all diagnostic tests, a definite clinical diagnosis should not be based on the results of a single test, but should only be made by the physician after all clinical and laboratory findings have been evaluated. The NSE test should not be used in cancer screening and should not replace any established clinical examination.
2. Studies have implicated possible interference in immunoassay results in some patients with known rheumatoid factor and antinuclear antibodies. Serum samples from patients who have received infusions containing mouse monoclonal antibodies for diagnostic or therapeutic purposes, may contain antibody to mouse protein (HAMA). Although we have added some agents to avoid the interferences, we cannot guarantee to eliminate all the effects of that.
3. Samples with NSE level above 200 ng/ml should be diluted to obtain accurate value.
4. High NSE values may be found in dialysis patients with leukaemic diseases.
5. Serum should not contain visible hemolysis since erythrocytes contain significant amounts of NSE.
6. Prolonged storage of whole blood can cause release of NSE from the blood cells.

### STORAGE

1. Unopened test kits should be stored at 2-8°C upon receipt and the microtiter plate should be kept in a sealed bag with desiccants to minimize exposure to damp air. The test kit may be used throughout the expiration date of the kit (one year from the date of manufacture). Refer to the package label for the expiration date.
2. Opened test kits will remain stable until the expiring date shown, provided it is stored as prescribed above.
3. A microtiter plate reader with a bandwidth of 10nm or less and an optical density range of 0-2.5 OD or greater at 450nm wavelength is acceptable for use in absorbance measurement.

### REFERENCES

1. Paus E, Nustad K. Immunoradiometric assay for  $\alpha\gamma$  and  $\gamma\gamma$ -Enolase (Neuron Specific Enolase), with use of Monoclonal antibodies and Magnetizable Polymer Particles. Clin.Chem, 35:2034, 1989.
2. Paus E., Risberg T., Establishment and Evaluation of A Radioimmunoassay for Neuron-Specific Enolase, Tumour Biol 10:23-30, 1989.
3. Cooper E.H., Pritchard J, bailey CC, Serum Neuron-Specific Enolase in children's cancer. Br J Cancer. 56:65-67,1987.
4. Wibe E., Paus E., Aamdal S., Neuron Specific Enolase (NSE) in serum of Patients with Malignant Melanoma, Cancer Letters, 52:29-31,1990.
5. Fossa S.D., Kiepp O., Paus E., Neuron-Specific (NSE) a serum tumor marker in seminoma. Br. J. Cancer, 65, 297-299,1992.
6. Skogseid I.M., nordby H.K., Urdal P., Paus E., Lilleaas F., Increased serum creatine kinase BB and Neuron-Specific Enolase following Head Injury indicates brain damage, Acta Neuronchir(Wien) 115:106-111,1992.
7. Pahlman S., Esscher T., Bergvall p. and Odelstad L. Purification and Characterization of human -specific Enolase, Radioimmunoassay Development Tumor Biol. 5,127-139, 1984.

### MANUFACTURER AND BRAND DETAILS

**ISO 13485:2016**

ISO 13485  
Quality  
Management for  
Medical Devices  
CERTIFIED

**Diagnostic Automation/Cortez Diagnostics, Inc.**  
21250 Califa Street, Suite 102 and 116,  
Woodland Hills, California 91367 USA

<b>Date Adopted</b>	2022-09
<b>Brand Name</b>	AccuDiag™
<b>REF</b> 6334-16	AccuDiag™ - NSE ELISA
<b>EC REP</b>	CEpartner4U, Esdoornlaan 13, 3951 DB Maarn, The Netherlands. <a href="http://www.cepartner4u.eu">www.cepartner4u.eu</a>
<b>Revision Date: 2016-01</b>	