



C-REACTIVE PROTEIN (CRP) LATEX TEST KIT

REF 310-100



WARNINGS AND PRECAUTIONS

1. Reagents containing sodium azide may combine with copper and lead plumbing to form highly explosive metal azides. Dispose of reagents by flushing with large amounts of water to prevent azide buildup.
2. For *in vitro* diagnostic use.
3. Positive and negative controls prepared using human sera found negative for hepatitis B surface antigen (HBsAg); however, handle controls as if potentially infectious.

REAGENT STORAGE AND PREPARATION

1. Reagents are stable until stated expiration date on bottle label when stored refrigerated (2 - 8°C).
2. **DO NOT FREEZE.**
3. The CRP latex reagent, once shaken must be uniform without visible clumping. When stored refrigerated slight sedimentation may occur and should be considered normal.
4. Do not use the latex reagent or controls if they become contaminated.

SPECIMEN COLLECTION AND STORAGE

1. Use fresh serum collected by centrifuging clotted blood.
2. If the test cannot be carried out on the same day, the serum may be stored between 2 - 8°C for no longer than 72 hours after collection.
3. For longer periods, the sample must be frozen.
4. As in all serological tests, hemolytic or contaminated serum must not be used.
5. Do not use plasma.

MATERIALS AND COMPONENTS

Materials provided with the test kit

1. CRP Latex Reagent
2. CRP Positive Control
3. CRP Negative Control
4. Glycine-Saline Buffer
5. Glass Slide
6. Stir Sticks

Materials required but not provided

1. Timer
2. Test Tubes
3. Test Tube Rack

ASSAY PROCEDURE

QUALITATIVE TEST:

1. Bring reagents and specimens to room temperature before use.
2. Place one drop (50 µl) of the CRP Positive Control on field #1 of the glass slide. Place one drop (50 µl) of the CRP Negative Control on field #2. The remaining fields are used for test specimens. Place one drop of the undiluted specimens on successive fields.
3. Gently resuspend the CRP Latex Reagent and add one drop to each test field. Use stir stick to spread reaction mixture over the entire test field.
4. Rotate the slide for three (3) minutes and read immediately under direct light.

INTENDED USE

C-Reactive protein (CRP) latex test is used for the qualitative and semi-quantitative measurement of C-reactive protein (CRP) in human serum.

SUMMARY AND EXPLANATION

C-Reactive Protein (CRP), the classic acute-phase of human serum, is synthesized by hepatocytes. Normally, it is present only in trace amounts in serum, but it can increase by as much as 1000-fold in response to injury or infection. The clinical measurement of CRP in serum, therefore, appears to be a valuable screening test for organic disease and a sensitive index of disease activity in inflammatory, infective, and ischemic conditions.^{1,2} MacLeod and Avery found that antibody produced against purified CRP provided a more sensitive test than the C-polysaccharide assay.³ Since that time a number of immunological assays have been devised to measure CRP such as capillary precipitation, double immunodiffusion, and radical immunodiffusion.^{4,5}

The CRP reagent kit is based on the principle of the latex agglutination assay described by Singer and Plotz.⁶ The major advantage of this method is the rapid three (3) minute reaction time.

ASSAY PRINCIPLE

The CRP reagent kit is based on an immunological reaction between CRP antisera bound to biologically inert latex particles and CRP in the test specimen. When serum containing greater than 0.8 mg/dl CRP is mixed with the latex reagent, visible agglutination occurs.

REAGENTS

1. **CRP Latex Reagent:** A suspension of uniform polystyrene particles coated with monospecific antihuman CRP (goat) in glycine buffer, pH 8.8 ± 0.5; reagent sensitivity adjusted to approximately 0.8 mg/dl. **MIX WELL BEFORE USING.**
2. **CRP Positive Control Serum:** A stabilized prediluted human serum containing more than 0.8 mg/dl CRP.
3. **CRP Negative Control Serum:** A stabilized prediluted human serum non-reactive with the test reagent.
4. **Glycine-Saline Buffer (20x):** pH 8.2 ± 0.1: A diluent containing 0.1M glycine and 0.15M NaCl. Dilute glycine-saline buffer before using by adding 1 part glycine saline buffer to 19 parts distilled water. All reagents contain 0.1% (w/v) sodium azide as preservative.



SEMI-QUANTITATIVE TEST:

1. Set up at least five test tubes: 1:2, 1:4, 1:8, 1:16, etc.
2. Dilute sample according to dilution factor on each test tube with diluted saline solution.
3. Place one drop of each of positive and negative controls onto separate slide fields. Place one drop of each dilution on successive fields of the reaction slides.
4. Gently resuspend the CRP Latex Reagent and add one drop to entire test field.
5. Mix well with the flat end of the stir stick. Gently rock the slide for three (3) minutes and read immediately under direct light.
6. The titer of the serum is the reciprocal of the highest dilution exhibiting a positive reaction multiplied by the concentration of the positive control.

mg/dL of serum =	concentration of positive control	x	reciprocal of last dilution showing a positive result
------------------	-----------------------------------	---	---

Dilution	Reciprocal	Conc. of serum (mg/dL)
1:2	2	1.6
1:4	4	3.2
1:8	8	6.4

3. It is important to determine the level of CRP for monitoring patient progress. This is due to (1) the concentration of CRP is an index of tissue damage incurred and (2) increasing or decreasing levels of CRP (e.g., daily) indicate the progress of inflammatory process.⁷

LIMITATIONS

1. Reaction time is critical. If reaction time exceeds three (3) minutes, drying of the reaction mixture may cause false positive results.
2. Freezing the CRP Latex Reagent will result in spontaneous agglutination.
3. Intensity of agglutination is not necessarily indicative of relative CRP concentration; therefore, screening reactions should not be graded.
4. A false negative can be attributed to a prozone phenomena (antigen excess). It is recommended, therefore, to check all negative sera by retesting at a 1:10 dilution with glycine buffer.

PERFORMANCE CHARACTERISTICS

1. Sensitivity: greater than 0.8 mg/dl
2. Linearity: 25 mg/dl.
3. Comparison:
 - A. Qualitative Results: A study performed using CRP Latex Reagent and a commercially available product yielded 100% accuracy.
 - B. Semi-quantitative Results: A panel of 32 CRP positive serum samples was assayed (semi-quantitated) on same day. The results of the study indicated that CRP Latex Reagent has 92.9 % precision.

QUALITY CONTROL

1. CRP Positive and Negative Control should be included in each test batch.
2. Acceptable performance is indicated when a uniform milky suspension with no agglutination is observed with the CRP Negative Control and agglutination with large aggregates is observed with the CRP Positive Control.

RESULTS

Negative Result: A negative reaction is indicated by a uniform milky suspension with no agglutination as observed with the CRP Negative Control.

Positive Result: A positive reaction is indicated by any observable agglutination in the reaction mixture. The specimen reaction should be compared to the CRP Negative Control (Figure 1).

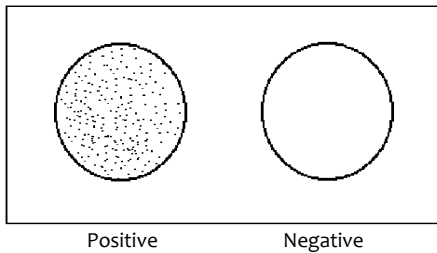


Figure 1.

REFERENCES

1. Pepys, M.B.. Lancet 1:653 (1981).
2. Werner, M.. Clin.Chem. Acta 25:299 (1969).
3. MacLeod, C.M., et. al.. J. Exp. Med 73:191 (1941).
4. Wood, HF., et. al.. J. Clin. Invest. 30: 616 (1951).
5. Mancini, G., et. al.. Immunochemistry 2:235 (1965).
6. Singer, J.M., et. al.. Am. J. Med 21: 888 (1956).
7. Fischer, C.L., Gill, C.W.. In Serum Protein Abnormalities. Boston, Little, Brown and Co., (1975).

Diagnostic Automation/Cortez Diagnostics, Inc. 21250 Califa Street, Suite 102 and 116, Woodland Hills, California 91367 USA	
Date Adopted	2022-09
REF 310-100	C-Reactive Protein (CRP) Latex test kit
EC REP	CEpartner4U, Esdoornlaan 13, 3951 DB Maarn, The Netherlands. www.cepartner4u.eu
Revision Date: 2018-12	

EXPECTED VALUES

1. CRP in healthy individuals is approximately 0.02-1.35mg/dl. The mean value in adults is 0.047mg/dl.
2. A weak positive correlation was found between CRP and age.