



DIAGNOSTIC AUTOMATION, INC.

23961 Craftsman Road, Suite E/F, Calabasas, CA 91302

Tel: (818) 591-3030 Fax: (818) 591-8383

onestep@rapidtest.com

technicalsupport@rapidtest.com

www.rapidtest.com

IVD



See external label



2°C-8°C



Σ=96 tests

REF

Cat # 4322Z

MICROWELL ELISA

FREE PROSTATE SPECIFIC ANTIGEN (f-PSA) ENZYME IMMUNOASSAY TEST KIT

Free PSA

Cat # 4322Z

Enzyme Immunoassay for the Quantitative Determination of Free Prostate Specific Antigen (f-PSA) in Human Serum

Intended Use

The f-PSA Enzyme Immunoassay test kit is intended for the quantitative determination of f-PSA in human serum.

Introduction

Human Prostate Specific Antigen (PSA) is a 33 kD serine proteinase which, in human serum, is predominantly bound to alpha 1-antichymotrypsin (PSA-ACT) and alpha 2-macroglobulin (PSA-AMG). Trace amounts of alpha 1-antitrypsin and inter-alpha trypsin inhibitor bound to PSA can also be found. Any remaining PSA is in the free form (f-PSA).^{1,3} Current methods of screening men for prostate cancer utilize the detection of the major PSA-ACT form. Levels of 4.0 ng/ml or higher are strong indicators of the possibility of prostatic cancer.⁴ However, elevated serum PSA levels have also been attributed to benign prostatic hyperplasia and prostatitis, leading to a large percentage of false positive screening results.⁵ A potential solution to this problem involves the determination of free PSA levels.^{1,3} Preliminary studies have suggested that the percentage of free PSA is lower in patients with prostate cancer than those with benign prostatic hyperplasia.⁵ Thus, the measurement of free serum PSA in conjunction with total PSA, can improve specificity of prostate cancer screening in selected men with elevated total serum PSA levels, which would subsequently reduce unnecessary prostate biopsies with minimal effects on cancer detection rates.⁶

Test Principle

The f-PSA EIA test is a solid phase two-site immunoassay. An anti-f-PSA monoclonal antibody is coated on the surface of the microtiter wells and another anti-PSA monoclonal antibody labeled with horseradish peroxidase is used as the tracer. The f-PSA molecules present in the standard solution or serum are "sandwiched" between the two antibodies. Following the formation of the coated antibody-antigen- antibody-enzyme complex, the unbound antibody-enzyme tracers are removed by washing. The horseradish peroxidase activity bound in the wells is then assayed by a colorimetric reaction. The intensity of the color formed is proportional to the concentration of f-PSA present in the sample.

Materials and Components

Materials provided with the test kit:

- Antibody-coated microtiter plate with 96 wells.

- Sample Diluent, 12 ml.
- Reference standards containing 0, 0.1, 0.5, 2.0, 5.0, and 10.0 ng/ml f-PSA, liquid, ready to use. 1 set.
- Enzyme Conjugate Reagent, 22 ml.
- TMB Substrate, 12ml
- Stop Solution, 12 ml.
- Wash Buffer Concentrate(50X),15ml.

Materials required but not provided:

- Precision pipettes: 0.10, 0.20, and 1.0 ml.
- Disposable pipette tips.
- Distilled water.
- Vortex mixer or equivalent.
- Absorbent paper or paper towels.
- Graph paper.
- A microtiter plate reader with a bandwidth of 10nm or less and an optical density range of 0-2 OD or greater at 450nm

Specimen Collection and Preparation

1. Blood should be drawn using standard venipuncture techniques and the serum should be separated from the red blood cells as soon as possible. Avoid grossly hemolytic, lipemic, or turbid samples.
2. Plasma samples collected in tubes containing EDTA, heparin, or oxalate may interfere with the test procedures and should be avoided.
3. Specimens should be capped and may be stored up to 48 hours at 2-8°C, prior to assaying. Specimens held for a longer time can be frozen at -20°C. Thawed samples must be mixed prior to testing.

Storage of test kits and instrumentation

1. Unopened test kits should be stored at 2-8°C upon receipt and the microtiter plate should be kept in a sealed bag with desiccants to minimize exposure to damp air. The test kit may be used throughout the expiration date of the kit (One year from the date of manufacture). Refer to the package label for the expiration date.
2. Opened test kits will remain stable until the expiring date shown, provided it is stored as prescribed above.
3. A microtiter plate reader with a bandwidth of 10nm or less and an optical density range of 0-2 OD or greater at 450nm wavelength is acceptable for use in absorbance measurement.

Reagent Preparation

1. All reagents should be brought to room temperature (18-22°C) and mixed by gently inverting or swirling prior to use. Do not induce foaming.
2. Dilute 1 volume of Wash Buffer (50x) with 49 volumes of distilled water. For example, Dilute 15 ml of Wash Buffer (50x) into distilled water to prepare 750 ml of washing buffer (1x). Mix well before use.

Assay Procedure

1. Secure the desired number of coated wells in the holder.
2. Dispense 100µl of standards, specimens, and controls into appropriate wells.
3. Dispense 100µl of sample diluent into each well.
4. Thoroughly mix for 10 seconds. It is very important to have a complete mixing in this step.
5. Incubate at 37°C for 60 minutes
6. Remove the incubation mixture by emptying plate contents into a suitable waste container.
7. Rinse and empty the microtiter wells 5 times with washing buffer (1X).
8. Strike the wells sharply onto absorbent paper or paper towels to remove all residual water droplets.
9. Dispense 200µl of Enzyme Conjugate Reagent into each well. Gently mix for 5 seconds.
10. Incubate at 37°C for 60 minutes.
11. Remove the incubation mixture by emptying plate contents into a suitable waste container.
12. Rinse and empty the microtiter wells 5 times with washing buffer (1X).
13. Strike the wells sharply onto absorbent paper to remove residual water droplets.
14. Dispense 100µl TMB solution into each well. Gently mix for 5 seconds.
15. Incubate at room temperature for 20 minutes in the dark.
16. Stop the reaction by adding 100µl of stop solution to each well.
17. Gently mix for 30 seconds to make sure that the blue color changes completely to yellow.

18. Using a microtiter plate reader, read the optical density at 450nm within 20 minutes.

Important Notes:

- The wash procedure is critical. Insufficient washing will result in poor precision and falsely elevated absorbance readings.
- It is recommended that if manual pipetting is used, no more than 32 wells be used for each assay run, since pipetting of all standards, specimens and controls should be completed within 3 minutes. A full plate of 96 wells may be used if automated pipetting is available.
- Duplication of all standards and specimens, although not required, is recommended.


Calculation of Results

Calculate the mean absorbance value (A_{450}) for each set of reference standards, controls, and patient samples. Construct a standard curve by plotting the mean absorbance obtained from each reference standard against its concentration in ng/ml on graph paper. The absorbance values are placed on the vertical, or Y-axis, and concentrations on the horizontal, or X-axis. Use the mean absorbance values for each specimen to determine the corresponding concentration of f-PSA in ng/ml from the standard curve.

Example of Standard Curve

Results of typical standard run with optical density reading at 450nm shown in the Y-axis against f-PSA concentrations shown in the X-axis. This standard curve is for the purpose of illustration only, and should not be used to calculate unknowns. Each user should obtain his or her own standard curve and data.

| f-PSA (ng/ml) | Absorbance (450nm) |
|---------------|--------------------|
| 0 | 0.006 |
| 0.1 | 0.032 |
| 0.5 | 0.155 |
| 2.0 | 0.597 |
| 5.0 | 1.302 |
| 10.0 | 2.361 |

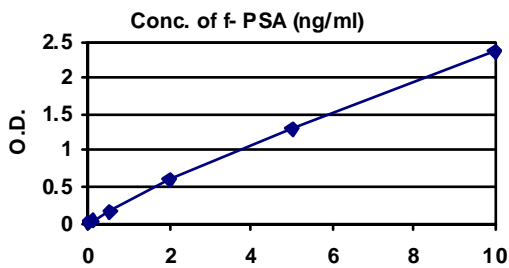

Sensitivity

The minimum detectable concentration of f-PSA in this assay is estimated to be 0.05 ng/ml.

Cross Reactivity

| Antigens | Concentration | % Cross-react. |
|---------------|---------------|----------------|
| PSA-ACT | 500 ng/ml | 0.2 |
| AFP | 10,000 ng/mL | 0 |
| CEA | 5,000 ng/mL | 0 |
| CA 125 | 1,000 U/mL | 0 |
| CA 15-3 | 1,000 U/mL | 0 |
| CA 19-9 | 1,000 U/mL | 0 |
| α -HCG | 1,000 ng/ml | 0 |
| β -HCG | 1,000 ng/mL | 0 |
| HCG | 50,000 mIU/ml | 0 |

Precision

| Intra-Assay | | | |
|-------------|------------|-------|------|
| | Replicates | S.D. | % CV |
| Level I | 20 | 0.005 | 13.1 |
| Level II | 20 | 0.011 | 4.4 |
| Level III | 20 | 0.128 | 3.2 |

Linearity

Two patient sera were serially diluted with 0 ng/ml standard. The average recovery was 108.7%.

| Sample A | | | |
|--------------------------|----------|----------|----------|
| Dilution | Expected | Observed | % Recov. |
| undiluted | 8.691 | 8.691 | 100 |
| 2X | 4.346 | 4.455 | 102.5 |
| 4X | 2.173 | 2.290 | 105.4 |
| 8X | 1.086 | 1.258 | 115.8 |
| 16X | 0.543 | 0.617 | 113.6 |
| 32X | 0.272 | 0.294 | 108.1 |
| 64X | 0.136 | 0.147 | 108.1 |
| Average Recovery: 107.6% | | | |

| Sample B | | | |
|--------------------------|----------|----------|----------|
| Dilution | Expected | Observed | % Recov. |
| undiluted | 7.015 | 7.015 | 100 |
| 2X | 3.508 | 3.516 | 100.2 |
| 4X | 1.754 | 1.834 | 104.6 |
| 8X | 0.877 | 0.970 | 110.6 |
| 16X | 0.438 | 0.509 | 116.2 |
| 32X | 0.219 | 0.249 | 113.7 |
| 64X | 0.110 | 0.136 | 123.6 |
| Average Recovery: 109.8% | | | |

Recovery

Equal parts of diluted patient sera were mixed to test for interference by unknown materials, such as drugs or hormones, in the assay. Concentrations of Free PSA were determined before (original and added) and after (observed). The average recovery was 99.2%.

| Sample 1 | | | | | |
|-------------------------|------------|-------|----------|----------|---------|
| Sampl e | Orig. Conc | Added | Expected | Observed | % Recov |
| A | 9.802 | 0.141 | 4.972 | 4.540 | 91.3 |
| B | 9.802 | 0.281 | 5.042 | 4.638 | 92.0 |
| C | 4.699 | 2.168 | 3.434 | 3.342 | 97.3 |
| D | 2.168 | 1.170 | 1.669 | 1.661 | 99.5 |
| E | 0.563 | 0.281 | 0.422 | 0.423 | 100.2 |
| F | 4.699 | 1.125 | 2.912 | 2.727 | 93.6 |
| G | 0.563 | 0.141 | 0.352 | 0.399 | 113.4 |
| Average Recovery: 98.2% | | | | | |

| Sample 2 | | | | | |
|--------------------------|------------|-------|----------|----------|---------|
| Sampl e | Orig. Conc | Added | Expected | Observed | % Recov |
| A | 6.825 | 0.098 | 3.462 | 3.301 | 95.3 |
| B | 6.825 | 0.195 | 3.510 | 3.156 | 89.9 |
| C | 3.342 | 1.563 | 2.453 | 2.503 | 102.0 |
| D | 1.706 | 0.781 | 1.244 | 1.205 | 96.9 |
| E | 0.391 | 0.195 | 0.293 | 0.327 | 111.6 |
| F | 3.342 | 0.781 | 2.062 | 1.964 | 95.2 |
| G | 0.391 | 0.098 | 0.245 | 0.271 | 110.6 |
| Average Recovery: 100.2% | | | | | |

Limitations of the Procedure

- Reliable and reproducible results will be obtained when the assay procedure is carried out with a complete understanding of the package insert instructions and with adherence to good laboratory practice.
- The wash procedure is critical. Insufficient washing will result in poor precision and falsely elevated absorbance readings.

3. Heterophilic antibodies such as human anti-mouse antibodies (HAMA) are frequently found in the serum of human subjects. Those antibodies can cause severe interference in many immunodiagnostic procedures. This assay has been designed to minimize that kinds of interference. Nevertheless, complete elimination of this interference from all patient specimens cannot be guaranteed. A test result that is inconsistent with the clinical picture and patient history should be interpreted with caution.

References

1. Christensson A, Bjork T, Nilsson O, et al. Serum Prostate Specific Antigen Complexed to α -1-Antichymotrypsin As An Indicator of Prostate Cancer. *J. of Urol.* 150:100-105; 1993.
2. Lilja H, Christensson A, Dahlen U, et al. Prostate-specific antigen in serum occurs predominantly in complex with α -1-antichymotrypsin. *Clin Chem.* 1991;37:1618-1625.
3. Stenman U-H, Leinonen J, Alfthan H, Rannikko S, Tuhkanen K, Alfthan O. A complex between prostate-specific antigen and α -1-antichymotrypsin is the major form of prostate-specific antigen in serum of patients with prostate cancer: assay of the complex improves clinical sensitivity for cancer. *Cancer Res.* 1991;51:222-226.
4. Catalona WJ, Smith DS, Ratliff TL, Basler JW. Detection of organ-confined prostate cancer is increased through prostate-specific antigen-based screening. *JAMA.* 1993;270:948-954.
5. Stamey TA, Yang N Hay AR, McNeal JE, Freiha, FS, Redwine E. Prostate-specific antigen as a serum marker for adenocarcinoma of the prostate. *N Engl J Med.* 1987;317:909-916.
6. Catalona et al. Percentage of Free Serum PSA and Prostate Cancer Detection. *JAMA.* 1995;274, No. 15: 1214-1220.
7. Catalona WJ, Smith DS, Ratliff TL, et al. Measurement of prostate-specific antigen in serum as a screening test for prostate cancer. *N Engl J Med.* 1991;324:1156-1161. Erratum: *N Engl J Med.* 1991;325:1324.
8. Carter HB, Pearson JD, Metter J, et al. Longitudinal evaluation of prostate-specific antigen levels in men with and without prostate disease. *JAMA.* 1992;267:2215-2220.
9. Smith DS, Catalona WJ. Rate of change in serum prostate-specific antigen levels as a method for prostate cancer detection. *J Urol.* 1994;152:1163-1167.
10. Benson MC, Whang IS, Pantuck A, et al. Prostate specific antigen density: a means of distinguishing benign prostatic hypertrophy and prostate cancer. *J Urol.* 1992; 147:815-816.
11. Catalona WJ, Hudson MA, Scardino PT, et al. Selection of optimal prostate specific antigen cut-offs for early detection of prostate cancer: receiver operating characteristic curves. *J Urol.* 1994;152:2037-2042.
12. Smith DS, Catalona WJ. The nature of prostate cancer detected through prostate specific antigen based screening. *J Urol.* 1994;152:1732-1736.

| Date Adopted | Reference No. |
|--------------|------------------|
| 2008-05-01 | DA-Free PSA-2008 |



DIAGNOSTIC AUTOMATION, INC.

23961 Craftsman Road, Suite E/F, Calabasas, CA
91302

Tel: (818) 591-3030 Fax: (818) 591-8383

ISO 13485-2003



Revision Date: 7/10/08

