

OneStep Strep A RapiCard™ InstaTest

REF 176550-1-44



A rapid test for the qualitative detection of Strep A antigens in throat swab specimens. For professional in vitro diagnostic use only.

INTENDED USE

The Cortez Strep A RapiCard™ InstaTest is a rapid chromatographic immunoassay for the qualitative detection of Strep A antigens from throat swab specimens to aid in the diagnosis of Group A Streptococcal infection.

SUMMARY AND EXPLANATION

Streptococcus pyogenes is non-motile gram-positive cocci, which contains the Lancefield group A antigens that can cause serious infections such as pharyngitis, respiratory infection, impetigo, endocarditis, meningitis, puerperal sepsis, and arthritis.¹ Left untreated, these infections can lead to serious complications, including rheumatic fever and peritonsillar abscess.² Traditional identification procedures for Group A Streptococci infection involve the isolation and identification of viable organisms using techniques that require 24 to 48 hours or longer.^{3,4}

The Strep A RapiCard™ InstaTest is a rapid test to qualitatively detect the presence of Strep A antigens in throat swab specimens, providing results within 5 minutes. The test utilizes antibodies specific for whole cell Lancefield Group A Streptococcus to selectively detect Strep A antigens in a throat swab specimen.

TEST PRINCIPLE

The Cortez Strep A RapiCard™ InstaTest a qualitative, lateral flow immunoassay for the detection of Strep A carbohydrate antigen in a throat swab. In this test, antibody specific to Strep A carbohydrate antigen is coated on the test line region of the test. During testing, the extracted throat swab specimen reacts with an antibody to Strep A that is coated onto particles. The mixture migrates up the membrane to react with the antibody to Strep A on the membrane and generate a color line in the test line region. The presence of this color line in the test line region indicates a positive result, while its absence indicates a negative result. To serve as a procedural control, a colored line will always appear in the control line region, indicating that proper volume of specimen has been added and membrane wicking has occurred.

REAGENTS

The test contains Strep A antibody coated particles and Strep A antibodies coated on the membrane.

MATERIALS AND COMPONENTS

Materials provided with the test kits

- Test Cassettes
- Workstation
- Extraction reagent 1 (2M NaNO₂)
- Positive control(Non-viable Strep A; 0.01% Proclin300)
- Negative control(Non-viable Strep C; 0.01% Proclin300)
- Extraction tubes
- Dropper tips
- Extraction reagent 2 (0.027M Citric acid)
- Sterile swabs
- Package insert

Materials required but not provided

- Timer

PRECAUTIONS

1. For professional in vitro diagnostic use only. Do not use after the expiration date.
2. Do not eat, drink or smoke in the area where the specimens and kits are handled.
3. Handle all specimens as if they contain infectious agents. Observe established precautions against microbiological hazards throughout the procedure and follow the standard procedures for proper disposal of specimens.
4. Wear protective clothing such as laboratory coats, disposable gloves and eye protection when specimens are assayed.
5. The used test should be discarded according to local regulations.
6. Humidity and temperature can adversely affect results.
7. Do not use test if pouch is damaged.
8. Reagent B contains an acidic solution. If the solution contacts the skin or eye, flush with large volumes of water.
9. The positive and negative controls contain Proclin300 as a preservative.
10. Do not interchange reagent bottle caps.
11. Do not interchange external control solution bottle caps.


SPECIMEN COLLECTION AND PREPARATION

1. Collect the throat swab specimen with the sterile swab that is provided in the kit. Transport swabs containing modified Stuart's or Amies medium can also be used with this product. Swab the posterior pharynx, tonsils and other inflamed areas. Avoid touching the tongue, cheeks and teeth with the swab.⁵
2. Testing should be performed immediately after the specimens have been collected. Swab specimens may be stored in a clean, dry plastic tube for up to 8 hours at room temperature or 72 hours at 2-8°C.
3. If a culture is desired, lightly roll the swab tip onto a Group A selective (GAS) blood agar plate before using the swab in the Cortez Strep A RapiCard™ InstaTest

ASSAY PROCEDURE

Allow the test, reagents, throat swab specimen, and/or controls to reach room temperature (15-30°C) prior to testing.

1. Remove the test cassette from the sealed foil pouch and use it within one hour. Best results will be obtained if the test is performed immediately after opening the foil pouch.
2. Hold the Extraction Reagent 1 bottle vertically and add 4 full drops (approximately 240 µL) of Extraction Reagent 1 to an extraction tube. Extraction Reagent 1 is red in color. Hold the Extraction Reagent 2 bottle vertically and add 4 full drops (approximately 160 µL) to the tube. Extraction Reagent 2 is colorless. Mix the solution by gently swirling the extraction tube. The addition of Extraction Reagent 2 to Extraction Reagent 1 changes the color of the solution from red to yellow. See illustration 1.
3. Immediately add the swab into the extraction tube, agitate the swab vigorously 15 times, Leave the swab in the extraction test tube for 1 minute. See illustration 2
4. Press the swab against the side of the tube and squeeze the bottom of the tube while removing the swab so that most of the liquid stays in the tube. Discard the swab. See illustration 3
5. Fit the dropper tip on top of the extraction tube. Place the test cassette on a clean and level surface. Add three drops of the solution (approx.100ul) to the sample well and then start the timer. Read the result at 5 minutes. Do not interpret the result after 10 minutes. See illustration 4 and illustration 5

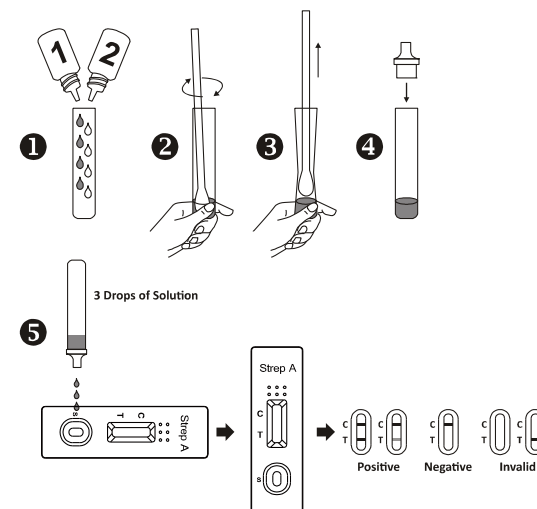

RESULTS

POSITIVE: Two lines appear. One colored line should be in the control line region (C) and another apparent colored line should be in the test line region (T). A positive result indicates that Strep A was detected in the specimen.

***NOTE:** The intensity of the color in the test line region (T) will vary depending on the concentration of Strep A present in the specimen. Therefore, any shade of color in the test line region (T) should be considered positive.




Diagnostic Automation / Cortez Diagnostics, Inc.

IMMUNO DIAGNOSTICS

NEGATIVE: One colored line appears in the control line region (C). No line appears in the test line region (T). A negative result indicates that Strep A antigen is not present in the specimen, or is present below the detectable level of the test. The patient's specimen should be cultured to confirm the absence of Strep A infection. If clinical symptoms are not consistent with results, obtain another specimen for culture.

INVALID: Control line fails to appear. Insufficient specimen volume or incorrect procedural techniques are the most likely reasons for control line failure. Review the procedure and repeat the test with a new test. If the problem persists, discontinue using the test kit immediately and contact your local distributor.

QUALITY CONTROL

Internal Quality Control

Internal procedural controls are included in the test. A colored line appearing in the control region (C) is an internal positive procedural control. It confirms sufficient specimen volume, adequate membrane wicking and correct procedural technique.

External Quality Control

It is recommended that a positive and negative external control be run every 25 tests, and as deemed necessary by internal laboratory procedures. External positive and negative controls are supplied in the kit. Alternatively, other Group A and non-Group A Streptococcus reference strains may be used as external controls. Some commercial controls may contain interfering preservatives; therefore, other commercial controls are not recommended.

Procedure for External Quality Control Testing

- Add 4 full drops of Extraction Reagent 1 and 4 full drops of Extraction Reagent 2 into an extraction tube. Tap the bottom of the tube gently to mix the liquid.
- Add 1 full drop of positive or negative control solution into the tube, holding the bottle upright.
- Place a clean swab into this extraction tube and agitate the swab in the solution by rotating it at least 15 times. Leave the swab in the extraction tube for 1 minute. Then express the liquid from the swab head by rolling the swab against the inside of the extraction tube and squeezing the extraction tube as the swab is withdrawn. Discard the swab.
- Continue with Step 5 of Directions For Use.
If the controls do not yield the expected results, do not use the test results. Repeat the test or contact your distributor.

PERFORMANCE CHARACTERISTICS

Sensitivity and Specificity

Using three medical centers for evaluation, a total of 526 throat swabs were collected from patients exhibiting symptoms of pharyngitis. Each swab was rolled onto a sheep blood agar plate, and then tested by the Strep A RapiCard™ InstaTest (Throat Swab). The plates were further streaked for isolation, and then incubated at 37 °C with 5-10% CO₂ and a Bacitracin disk for 18-24 hours. The negative culture plates were incubated for an additional 18-24 hours. Possible GAS colonies were subcultured and confirmed with a commercially available latex agglutination grouping kit. Of the 525 total specimens, 401 were confirmed to be negative and 124 were confirmed to be positive by culture. During this study, one Strep F specimens yielded positive results with the Test. One of these specimens was re-cultured, then re-tested and yielded a negative result. Three additional different Strep F strains were cultured and tested for cross-reactivity and also yielded negative results.

Method	Results	Culture		Total Results
		Positive	Negative	
Strep A RapiCard™				

InstaTest	Positive	116	9	125
	Negative	6	395	401
Total Results		122	404	526

Relative Sensitivity: 95.1% (95%CI*: 89.6%-98.2%) *Confidence Interval
Relative Specificity: 97.8% (95%CI*: 95.8%-99%)
Accuracy: 97.1% (95%CI*: 95.3%-98.4%)

Positive Culture Classification	Strep A Rapid Test/Culture	% Agreement
Rare	8/10	80.0%
1+	18/20	90.0%
2+	19/20	95.0%
3+	33/34	97.1%
4+	38/38	100.0%

Cross Reactivity

The following organisms were tested at 1.0 x 10⁷ organisms per test and were all found to be negative when tested with the Strep A RapiCard™ InstaTest. No mucoid-producing strains were tested.

Group B Streptococcus	Neisseria meningitidis	Serratia marcescens
Group F Streptococcus	Neisseria sicca	Klebsiella pneumoniae
Streptococcus pneumoniae	Branhamella catarrhalis	Bordetella pertussis
Streptococcus mutans	Group C Streptococcus	Neisseria gonorrhoea
Staphylococcus aureus	Group G Streptococcus	Neisseria subflava
Corynebacterium diphtheria	Streptococcus sanguis	Hemophilus influenza
Candida albicans	Staphylococcus epidermidis	Pseudomonas aeruginosa
Enterococcus faecalis		

LIMITATIONS OF PROCEDURE

- The Cortez Strep A RapiCard™ InstaTest is for in vitro diagnostic use only. The test should be used for the detection of Strep A antigen in throat swab specimens only. Neither the quantitative value nor the rate of increase in Strep A antigen concentration can be determined by this qualitative test.
- This test will only indicate the presence of Strep A antigen in the specimen from both viable and non-viable Group A Streptococcus bacteria.
- A negative result should be confirmed by culture. A negative result may be obtained if the concentration of the Strep A antigen present in the throat swab is not adequate or is below the detectable level of the test.
- Excess blood or mucus on the swab specimen may interfere with test performance and may yield a false positive result. Avoid touching the tongue, cheeks, and teeth⁵ and any bleeding areas of the mouth with the swab when collecting specimens.
- As with all diagnostic tests, all results must be interpreted together with other clinical information available to the physician.

EXPECTED VALUES

Approximately 15% of pharyngitis in children ages 3 months to 5 years is caused by Group A beta hemolytic Streptococcus.⁶ In school-aged children and adults, the incidence of Strep throat infection is about 40%.⁷ This disease usually occurs in the winter and early spring in temperate climates.³

STORAGE

Store as packaged in the sealed pouch at room temperature or refrigerated (2-30°C). The test is stable through the expiration date printed on the sealed pouch. The test must remain in the sealed pouch until use. DO NOT FREEZE. Do not use beyond the expiration date.



REFERENCES

- Murray, P.R., et al. Manual of Clinical Microbiology, 6th Edition, ASM Press, Washington D.C., 1995, p. 299-307.
- Webb, KH. Does Culture Confirmation of High-sensitivity Rapid Streptococcal Tests Make Sense? A Medical Decision Analysis. Pediatrics (Feb 1998), 101:2, 2.
- Bisno AL, Gerber MA, Gwaltney JM, Kaplan EL, Schwartz RH. Diagnosis and Management of Group A Streptococcal Pharyngitis. Clinical Infectious Diseases (1997), 25: 574-83.
- Needham CA, McPherson KA, Webb KH. Streptococcal Pharyngitis: Impact of a High-sensitivity Antigen Test on Physician Outcome. Journal of Clinical Microbiology (Dec 1998), 36: 3468-3473.
- Shea, Y.R., Specimen Collection and Transport, Clinical Microbiology Procedures Handbook, Isenberg, H.D., American Society of Microbiology, Washington D.C., 1.1.1-1.1.30, 1992.
- Nussinovitch, M, Finkelstein Y, Amir J, Varsano, I. Group A beta-hemolytic streptococcal pharyngitis in preschool children aged 3 months to 5 years. Clinical Pediatrics (June 1999), 38: 357-360.
- Woods WA, Carter CT, Stack M, Connors Jr AF, Schlager TA. Group A Streptococcal Pharyngitis in Adults 30 to 65 years of age. Southern Medical Journal (May 1999), 491-492.

Diagnostic Automation / Cortez Diagnostics, Inc.
21250 Califa St, Suite 102 and 116, Woodland Hills, CA 91367 USA Phone: 818-591-3030, Fax: 818-591-8383
Email: onestep@rapidtest.com Website: www.rapidtest.com



Diagnostic Automation / Cortez Diagnostics, Inc.
I M M U N O D I A G N O S T I C S

<p>ISO 13485 ISO 9001</p> 	
 <p>Diagnostic Automation/ Cortez Diagnostics, Inc. 21250 Califa St, Suite 102 and 116, Woodland Hills, California 91367 USA</p>	
Date Adopted	2017-03-30
REF 176550-1-44	CORTEZ - OneStep Strep A RapiCard™ InstaTest
EC REP	<p>CEpartner4U, Esdoornlaan 13, 3951DB Maarn. The Netherlands. www.cepartner4u.eu</p>
Revision Date: 2016-11-16	

Diagnostic Automation/ Cortez Diagnostics, Inc.
 21250 Califa St, Suite 102 and 116, Woodland Hills, CA 91367 USA Phone: 818-591-3030, Fax : 818-591-8383
 Email: onestep@rapidtest.com Website: www.rapidtest.com