

**OneStep
 Tramadol Urine
 RapiDip™ InstaTest**

REF 121118-1



Sensitivity 200 ng/ml

INTENDED USE

The Cortez OneStep Tramadol RapiDip™ InstaTest are immunochromatography based one step in vitro test.

SUMMARY AND EXPLANATION

Tramadol is a quasi-narcotic analgesic used in the treatment of moderate to severe pain. It is a synthetic analog of codeine, but has a low binding affinity to the mu-opioid receptors. Large doses of tramadol can develop tolerance and physiological dependency and lead to its abuse. Tramadol is extensively metabolized after oral administration. Approximately 30% of the dose is excreted in the urine as unchanged drug, whereas 60% is excreted as metabolites. The major pathways appear to be N- and O- demethylation, glucuronidation or sulfation in the liver.

TEST PRINCIPLE

The Cortez OneStep Tramadol RapiDip™ InstaTest device is based on the principle of specific immunochemical reaction between antibodies and antigen to analyze particular compound in human urine specimen. The assay relies on the competition for binding antibody. When drug is present in the urine specimen, it competes with drug conjugate for the limited amount of antibody-dye conjugate. When the amount of drug is equal or more than the cut-off, it will prevent the binding of drug conjugate to the antibody. Therefore, a positive urine specimen will not show a colored band

on the test line zone, indicating a positive result, while the presence of a colored band indicates a negative result.

A control line is present in the test window to work as procedural control. This colored band should always appear on the control line zone if the test device is stored in good condition and the test is performed appropriately.

SPECIMEN COLLECTION AND PREPARATION

Fresh urine does not require any special handling or pretreatment. Specimen should be collected in a clean, dry, plastic or glass container. If the assay is not performed immediately, urine specimen may be refrigerated at 2-8 °C or frozen up to 7 days. Specimens should be brought to room temperature before testing. Urine specimens exhibiting a large amount of precipitate or turbidity should be centrifuged or allowed to settle before testing. Avoid contact with skin by wearing gloves and proper laboratory attire.

MATERIALS AND COMPONENTS

Materials provided with the test kits

1. The Cortez OneStep Tramadol RapiDip™ InstaTest device. The amount of each coated antigen and/or antibody on the strip is less than 1.0 mg for antigen conjugate and is less than 1.0 mg for goat anti-mouse IgG antibody.
2. Test zone: contains drug bovine protein antigen conjugates.
3. Control zone: contains Goat anti-mouse IgG antibody.
4. Conjugate pad: contains mice monoclonal anti-drug antibody.
5. Instruction for use.

Materials required but not provided

1. Urine collection container.
2. Timer or clock.

ASSAY PROCEDURE


1. Bring all materials and specimens to room temperature.
2. Remove the test strip from the sealed foil pouch.
3. Dip the strip into the urine specimen with the arrow pointing toward the sample. The sample level should not be higher than the arrow pointed maximum line.

4. Hold the strip in the urine until a reddish color appears at the test area (approximately 20 seconds).
5. Withdraw the strip and place it face up on a clean, non-absorptive surface or leave the strip in urine if the urine level is not higher than arrow pointed maximum line.
6. Read the results at 5 minutes after adding the sample.

RESULTS

- **Negative:** Two colored bands form. The appearance of two colored bands, one in test line zone and the other in control line zone, indicates negative result for that particular test(s). The negative result does not indicate the absence of drug in the specimen. It only indicates the level of tested drug in the specimen is less than cut-off level.
- **Positive:** One colored band forms. One colored band appears in control line zone. No colored band is found in test line zone. This is an indication that the level of tested drug(s) in the specimen is above the cut-off level.
- **Invalid:** If there are no colored bands in control line zone, the test result is invalid. Retest the sample with a new device.

Note: A borderline (+/-) in test line zone be considered negative result.

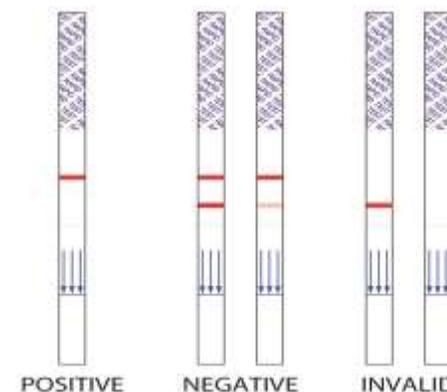

PERFORMANCE CHARACTERISTICS

A Accuracy

The accuracy of the The Cortez OneStep Tramadol RapiDip™ InstaTest panels were evaluated in each component strip and in comparison to GC/MS method at the following concentration: Tramadol 200 ng/ml (TRA). The results of each component strip are listed below:

Tramadol The accuracy of the Tramadol test was evaluated in comparison to GC/MS at a cut-off of 200 ng/ml of tramadol. Eighty one (81) urine specimens with GC/MS confirmed tramadol concentration were evaluated in this study. The results are summarized and presented below:
 Positive % agreement: 95, Negative % agreement: 98

B Sensitivity

The cut-off concentrations (sensitivity level) of The Cortez OneStep Tramadol RapiDip™ InstaTest are determined to be: TRA 200 ng/ml.

C Precision

The precision of The Cortez OneStep Tramadol RapiDip™ InstaTest were determined by conducting the test with spiked controls and interpreted the results by three individuals to verify the random error of visual interpretation. The results of 40 samples each of 50% above and 50% below cut-off specimens are 100% agreed by three observers. The test results were found to have no significant differences between these three observers.

D Specificity

The specificity for The Cortez OneStep Tramadol RapiDip™ InstaTest were tested by adding various drugs, drug metabolites, and other compounds that are likely to be present in urine. All compounds were prepared in drug-free normal human urine.

1. Interference Testing

The The Cortez OneStep Tramadol RapiDip™ InstaTest performance at cut-off level is not affected when pH and Specific Gravity ranges of urine specimen are at 4.0 to 9.0 and 1.005 to 1.035

The following substances were tested and confirmed did not interfere with The Cortez OneStep Tramadol RapiDip™ InstaTest at the listed concentrations.

Glucose	2000 mg/dl,
Human albumin	2000 mg/dl
Human hemoglobin	10 mg/dl,
Urea	4000 mg/dl
Uric acid	10 mg/dl

2. Specificity

The following table lists compounds that are detected by Cortez OneStep Tramadol RapiDip™ InstaTest which produced positive results when tested at levels equal or greater than the concentrations listed below:

Tests	Compounds	Cut-off (ng/ml)
Tramadol	Tramadol	200
	N-desmethyl-tramadol	500
	O-desmethyl-tramadol	20,000

The following compounds show no cross-reactivity at concentration up to 100 µg/ml unless specified.

Acetaminophen	4-Acetamidophenol	Acetylsalicylic acid
Arterenol	Aspartame	Ascorbic acid
Caffeine	Camphor	Chloroquine
Cortisone	Deoxyephedrine	Dextromethorphan
Digoxin	Diphenhydramine	Ecgonine
Ephedrine	Epinephrine	Gentisic acid
Histamine	Hydrochlorothiazide	Homatrophine
Isoproterenol	Lidocaine	Meperidine
Methylphenidate	Neomycin	Niacinamide
Penicillin G	Phenylethylamine-α	Phenylpropanolamine
Pseudoephedrine	Quinine antidine	Salicylic acid
Tetrahydrozoline	Theophylline	Thioridazine
Tryptophan	Tyramine	

Amikacin
 Atrophine
 Chlorpheniramine
 Digitoxin
 Ecgonine methyl ester
 Guaiacol glycer ester
 Ibuprofen
 Methaqualon
 Perphenazine
 Promethazine
 Tetracycline
 Trifluoperazine

Good Laboratory Practice recommends the daily use of control materials to validate the reliability of device. Control materials should be assayed as clinical specimen and challenging to the assay cutoff concentration, e.g., 25% above and below cutoff concentration. If control values do not fall within establish range, assay results are invalid. Control materials which are not provided with this test kit are commercially available.

The Cortez OneStep Tramadol RapiDip™ InstaTest provides a built-in process control with a different antigen/antibody reaction at the control region (C). This control line should always appear regardless the presence of drug or metabolite. If the control line

does not appear, the test device should be discarded and the obtained result is invalid. The presence of this control band in the control region serve as 1) verification that sufficient volume is added, 2) that proper flow is obtained.

LIMITATION OF PROCEDURE

The assay is designed for use with human urine only. A positive result with any of the tests indicates only the presence of a drug/metabolite and does not indicate or measure intoxication. There is a possibility that technical or procedural error as well other substances in certain foods and medicines may interfere with the test and cause false results. Please refer "SPECIFICITY" section for lists of substances that will produce either positive results, or that do not interfere with test performance. If a drug/metabolite is found present in the urine specimen, the assay does not indicate frequency of drug use or distinguish between drug of abuse and certain foods and medicines.

EXPECTED VALUES

The Cortez OneStep Tramadol RapiDip™ InstaTest is a qualitative assay. It identifies the drug(s) in human urine at its cut-off concentration or higher. The concentration of the drug(s) cannot be determined by this assay. The test is intended to distinguish negative result from presumptive positive result. All positive results must be confirmed using an alternate method, preferably GC/MS.

PRECAUTION

1. For in vitro diagnostic and forensic use only.
2. Do not use the product beyond the expiration date.
3. Handle all specimens as potentially infectious.
4. Humidity sensitive product, do not open foil pouch until it is ready to be tested.
5. Use a new urine specimen cup for each sample to avoid cross contamination.
6. The test device should be stored at 2 to 30°C and will be effective until the expiration date stated on the package. The product is humidity-sensitive and should be used immediately after being open. Any improperly sealed product should be discarded.

ISO 13485 ISO 9001 	
 Diagnostic Automation/ Cortez Diagnostics, Inc. 21250 Califa St, Suite 102 and 116, Woodland Hills, California 91367 USA	
Date Adopted	2016-01-28
 121118-1	CORTEZ- OneStepm Tramadol Urine RapiDip™ InstaTest
Revision Date: 2007-08-27	