

Diagnostic Automation/Cortez Diagnostics, Inc.

AccuDiag[™] CA-15-3 **ELISA Kit**

REF 6333-16

IVD See External Label 2°C

CA-15-3 ELISA	
Method	Enzyme Linked Immunosorbent Assay
Principle	Sandwich Complex
Detection Range	0-240U/mL
Sample	100 μL Serum
Specificity	97.5%
Sensitivity	5 U/mL
Incubation Time	140 minutes
Shelf Life	12 Months from the manufacturing date

PRODUCT FEATURES



Very easy to use with little training



Highly specific and consistent Assay



Provides accurate results quickly



Reading of results both visually and as absorbance data

Enzyme Immunoassay for the Quantitative Measurement of Breast Cancer Antigen (CA15-3) in Human Serum.

INTENDED USE

The CA-15-3 ELISA test is intended for use as a monitoring and screening test for breast cancer. An abnormal result (i.e., elevated serum CA-15-3 level) indicates further clinical management. CA-15-3 is a useful tumor marker for patients in clinical remission following treatment. Post-operative serum CA-15-3 values which fail to return to normal, strongly suggest the presence of residual tumor, while tumor recurrence is often accompanied by a rise of serum CA-15-3 before progressive disease is clinically evident.

SIGNIFICANCE AND SUMMARY

Breast cancer is the most common life-threatening malignant lesion in women of many developed countries today, with approximately 180,000 new cases diagnosed every year. Roughly half of these newly diagnosed patients are node-negative; however, 30% of these cases progress to metastatic disease. There are a number of tumor markers that can help clinicians to identify and diagnose which breast cancer patients will have an aggressive disease and which will have an indolent course. These markers include estrogen and progesterone receptors, DNA ploidy and percent-S phase profile, epidermal growth factor receptor, HER-2/neu oncogene, p53 tumor suppressor gene, cathepsin D, proliferation markers and CA-15-3. CA-15-3 is most useful for monitoring patients post-operatively for recurrence, particularly metastatic diseases. 96% of patients with local and systemic recurrence have elevated CA-15-3, which can be used to predict recurrence earlier than radiological and clinical criteria. A 25% increase in the serum CA-15-3 is associated with progression of carcinoma. A 50% decrease in serum CA-15-3 is associated with response to treatment. CA-15-3 are more sensitive than CEA in early detection of breast cancer recurrence. In combination with CA-12-5, CA-15-3 has been shown to be useful in early detection of relapse of ovarian cancer. CA-15-3 levels are also increased in colon, lung and hepatic tumors.

ASSAY PRINCIPLE

The CA-15-3 ELISA test is a two-site solid phase enzyme immunoassay. The molecules of CA-15-3 are "sandwiched" between two monoclonal antibodies. One coated to the bottom of the wells of the microtiter plates and the other linked to the horseradish peroxidase (enzyme conjugate). After incubation and washing, the enzymatic reaction develops a color which is proportional to the amount of CA-15-3 molecules present in the assay.

SPECIMEN COLLECTION AND PREPARATION

- Blood should be drawn using standard venipuncture techniques and the serum should be separated from the red blood cells as soon as practical. Avoid grossly hemolytic, lipemic or turbid samples.
- Plasma samples collected in tubes containing EDTA, heparin, or oxalate may interfere with test procedures and should be avoided. Specimens should be capped and may be stored for 48 hour at 2-8°C prior to assaying. Specimens held for a longer time can be frozen at -20°C for 6 months prior to assaying. Thawed samples should be inverted several times to mix prior to testing.

MATERIALS AND COMPONENTS

Materials provided with the test kit

- Monoclonal Anti-CA-15-3 antibody coated microtiter plate with 96 wells. 1.
- 2. Sample diluent, 100 ml.
- Enzyme conjugate reagent, 12 ml. 3.
- CA-15-3 reference standard set, containing o, 15, 30, 60, 120, and 240 Unit/ml (liquid, ready for use) or lyophilized form.
- Wash Buffer Concentrate (50X), 15 ml.
- 6. TMB Substrate, 12 ml.
- Stop solution, 12 ml. 7.
- Control set (optional)

Materials required but not provided

- Precision pipettes and tips, 0.1 ml, 0.2 ml, 1 ml, and 5 ml.
- Distilled water.
- Disposable pipette tips. 3.

Diagnostic Automation/Cortez Diagnostics, Inc.

21250 Califa St, Suite 102 and 116, Woodland Hills, CA 91367 USA Phone: 818-591-3030, Fax: 818-591-8383

Email: onestep@rapidtest.com Website: www.rapidtest.com

6333-P1



Diagnostic Automation/Cortez Diagnostics, Inc.

((

- Vortex mixer.
- 5. Absorbent paper or paper towel.
- 6. A microtiter plate reader at 450nm wavelength, with a bandwidth of 10nm or less and an optical density range of 0-2-5 OD or greater.
- 7. Graph paper.

REAGENT PREPARATION

- All reagents should be brought to room temperature (18-22°C) before use. All reagents should be mixed by gently inverting or swirling prior to use. Do not induce foaming.
- If reference standards are lyophilized, reconstitute each standard with 0.5 ml distilled water. Allow the reconstituted material to stand for at least 20 minutes. Reconstituted standards should be sealed and stored at 2-8°C.
- 3. Dilute 1 volume of Wash Buffer Concentrate (50x) with 49 volumes of distilled water. For example, dilute 15 ml of Wash Buffer Concentrate (50x) into 735 ml of distilled water to prepare 750 ml of washing buffer (1x). Mix well before use.

ASSAY PROCEDURE

Important Note:

The CA-15-3 standards have already been prediluted and are ready for use. Please DO NOT dilute again!

- 1. Patient serum and control serum should be diluted, 51-fold, before use. Prepare a series of small tubes (such as 1.5 ml microcentrifuge tubes) and mix 20 µl serum with 1.0 ml Sample Diluent.
- 2. Secure the desired number of coated wells in the holder. Dispense 100 µl of CA-15-3 standards, diluted specimens, and diluted controls into the appropriate wells. Gently mix for 10 seconds.
- 3. Incubate at 37°C for 1 hour.
- 4. Remove the incubation mixture by emptying the plate content into a waste container. Rinse and empty the microtiter plate 5 times with 1x wash buffer. Strike the microtiter plate sharply onto absorbent paper or paper towels to remove all residual water droplets.
- 4. Dispense 100 μ l of enzyme conjugate reagent into each well. Gently mix for 10 seconds
- 5. Incubate at 37°C for 1 hour.
- 6. Remove the contents and wash the plate as described in step 4 above.
- Dispense 100 μl TMB substrate reagent into each well. Gently mix for 10 seconds. Incubate at room temperature for 20 minutes. Stop the reaction by adding 100 μl of Stop Solution to each well.
- Gently mix for 10 seconds ensuring that the blue color completely changes to yellow. Read the optical density at 450nm with a microtiter plate reader within 15 minutes.

Important Note:

- The wash procedure is critical. Insufficient washing will result in poor precision and falsely elevated absorbance readings.
- It is recommended that no more than 32 wells be used for each assay run,
 if manual pipetting is used, since pipetting of all standards, specimens
 and controls should be completed within 5 minutes. A full plate of 96
 wells may be used if automated pipetting is available.
- Duplication of all standards and specimens, although not required, is recommended.


RESULTS

Calculate the mean absorbance value for each set of CA-15-3 reference standards, specimens and controls. Construct a standard curve by plotting the mean absorbance obtained from each reference standard against its

concentration in units per ml on linear graph paper, with absorbance values on the vertical or Y axis and concentrations on the horizontal or X axis. Use the mean absorbance values for each specimen to determine the corresponding concentration of CA-15-3 in units per ml from the standard curve. It is recommended that samples be analyzed in duplicates.

Since the CA-15-3 standards have already been diluted 51-fold, there is no need for the samples or controls to be multiplied by the dilution factor.

CA-15-3 Values (U/ml)	Absorbance (450 nm)
0	0.021
15	0.425
30	0.693
60	1.214
120	1.956
240	2.845

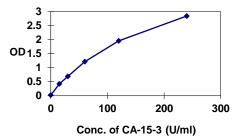

EXPECTED VALUES AND SENSITIVITY

Healthy women are expected to have CA-15-3 values below 35 U/ml. The minimum detectable concentration of CA-15-3 in this assay is estimated to be 5 U/ml.

LIMITATIONS OF THE PROCEDURE

There are some limitations of the assay:

- As with all diagnostic tests, a definite clinical diagnosis should not be based on the results of a single test, but should only be made by the physician after all clinical and laboratory findings have been evaluated.
- 2. Studies have implicated possible interference in immunoassay results in some patients with known rheumatoid factor and antinuclear antibodies. Serum samples from patients who have received infusions containing mouse monoclonal antibodies for diagnostic or therapeutic purposes, may contain antibody to mouse protein (HAMA). Although we have added some agents to avoid the interferences, we cannot guarantee it will eliminate all the effects of that.

STORAGE

- Unopened test kits should be stored at 2-8°C upon receipt and the microtiter plate should be kept in a sealed bag with desiccants to minimize exposure to damp air. The test kit may be used throughout the expiration date of the kit (one year from the date of manufacture). Refer to the package label for the expiration date.
- Opened test kits will remain stable until the expiring date shown, provided it is stored as prescribed above.

Diagnostic Automation/Cortez Diagnostics, Inc.

21250 Califa St, Suite 102 and 116, Woodland Hills, CA 91367 USA Phone: 818-591-3030, Fax: 818-591-8383

Email: onestep@rapidtest.com Website: www.rapidtest.com

6333-P1 Page 2 of 3



\ <u>\</u>

Diagnostic Automation/Cortez Diagnostics, Inc.


GNOSTICS

 A microtiter plate reader with a bandwidth of 10nm or less and an optical density range of 0-2.5 OD or greater at 450nm wavelength is acceptable for use in absorbance measurement.

REFERENCES

- Aziz DC, Rittenhouse HJ, Ranken R. Use and interpretation of tests in oncology. Santa Monica: Specialty Laboratories, 1991.
- Aziz DC. Quantitation of estrogen and progesterone receptors by immunocytochemical and image analyses. A J Clin Pathol 1992;98:105-11
- Aziz DC, Peter JB. DNA ploidy and cell-cycle analysis. Tools for assessment of cancer prognosis. J Clin Pathol 1991;5:422-38.
- Clark GM, Dressler LG, Owens MA, Dounds G, Oldaker T, McGuire WL. Prediction of relapse or survival in patients with node-negative breast cancer by DNA flow cytometry. N Engl J Med 1989;320:627-33.
- Elledge RM, McGuire WL. Prognostic factors and therapeutic decisions in axillary node-negative breast cancer. Annu Rev Med 1993;44:201-10.
- 6. Foekens JA, Rio C, Seguin P, et al. Prediction of relapse and survival in breast cancer patients by pS2 protein. **Cancer Res** 1990; 50-3832-7.
- Isola J, Visakorp T, Holli K, Kallionieml D. Association of p53 expression with other prognostic factors and long-term survival in node-negative breast cancer. J Cell Biochem 1992;(Suppl 16D):101.
- 8. Kute TE, Shao ZM, Snugg NK, Long RT, Russell GB, Case LD. Cathepsin D as a prognostic indicator for node-negative breast cancer patients using both immunoassays and enzymatic assays. **Cancer Res** 1992;52-198-203.
- McGuire WL, Tandon AK, Allred D, Chamnes GC, Clark GM. How to use prognostic factors in axillary node a) negative breast cancer patients. J Natl Cancer Inst 1990;82:1006-7.
- Nicholson S, Richard J, Sainsbury C, et al. Epidermal growth factor receptor (EGFr): results of a 6 year follow up study in operable breast cancer with emphasis on the node-negative subgroup. Br J Cancer 1991;63:146-50.
- Somerville JE, Clarke LA, Biggart JD. C-erb B-2 overexpression and histological type of in-situ and invasive breast carcinoma. J Clin Pathol 1992;45-16-20.
- 12. Ueronese S, Gambacorta M. Detection of Ki-67 rate in breast cancer. **Am J Clin Patho**l 1991;95:30-4.
- 13. Lotnicker M, Pavesi F, Scarabelli M. Tumor associated antigens CA15-3 and CA-125 in ovarian cancer. Int. J. Biolog Markers 1991; 6:115

MANUFACTURER AND BRAND DETAILS

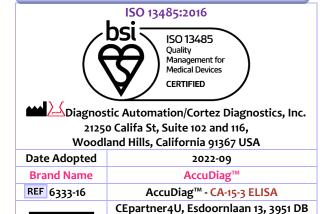

Revision Date: 2016-01

Maarn, The Netherlands. www.cepartner4u.eu

Diagnostic Automation/Cortez Diagnostics, Inc.

21250 Califa St, Suite 102 and 116, Woodland Hills, CA 91367 USA Phone: 818-591-3030, Fax: 818-591-8383

Email: onestep@rapidtest.com Website: www.rapidtest.com

EC

REP