

**AccuDiag™
 Zearalenone
 ELISA Kit**

REF 5148-8



Test	Zearalenone ELISA
Recovery (spiked samples)	86 - 110 %
Incubation Time	45 min.
Shelf Life	12 Months from the manufacturing date
Sensitivity	0.04 ng/mL

INTENDED USE

Zearalenone in addition to fumonisin, deoxynivalenol and other trichothecenes belongs to the fusarium toxins. These toxins are already produced on the field in consequence of a contact of the cereals by fusarium species. These moulds infect grain and other types of food like peanuts and beans already during their growth.

When a considerable amount of zearalenone contaminated feed is taken up by cows, it can also be detected in their milk. Even in beer it could be found. Zearalenone shows a strong estrogen-like activity. Thus zearalenone can cause an enlargement of the uterus, diminution of the ovarian glands and even infertility. Zearalenone is one of the main contaminants of farm products, which can be taken up by humans and animals.

In the European Union the limits are 20 – 400 ppb for food products. Thus a monitoring of food and feed with respect to the concentration of zearalenone is obligatory.

The Diagnostic Automation Inc. (DAI) Zearalenone ELISA represents a very highly sensitive detection system and is particularly capable of the quantification of zearalenone contaminations in cereals, milk, nuts, fruit and beer/gyle.

TEST PRINCIPLE

The Diagnostic Automation, Inc. quantitative test is based on the principle of the enzyme-linked immunosorbent assay. An antibody directed against zearalenone is coated on the surface of a microtiter plate. Zearalenone containing samples or standards and a zearalenone-peroxidase conjugate are given into the wells of the microtiter plate. The conjugate competes with the zearalenone of the samples/standards for the limited number of antibody sites. After 30 minutes incubation at room temperature the wells are washed with diluted washing solution to remove unbound material. A substrate solution is added and incubated for 15 minutes, resulting in the development of a blue color. The color development is inhibited by the addition of a stop solution, and the color turns yellow. The yellow color is measured photometrically at 450 nm. The concentration of zearalenone is inversely proportional to the color intensity of the test sample.

MATERIALS AND COMPONENTS

Materials provided with the test kits

The kit contains reagents for 96 determinations. They have to be stored at 2-8°C. Expiry data are found on the labels of the bottles and the outer package.

1. **Microtiter plate** consisting of 12 strips with 8 breakable wells each, coated with anti-zearalenone antibody.
2. **Zearalenone Standards (0; 0.05; 0.15; 0.5; 1.5; 3.5 ng/mL)**: 6 vials with 1.0 mL each, dyed red, ready-to-use.
3. **Conjugate (Zearalenone-peroxidase)**: 6 mL, dyed red, ready-to-use.
4. **Substrate Solution (TMB)**: 15 mL; ready-to-use.
5. **Stop Solution (0.5 M H₂SO₄)**: 15 mL; ready-to-use.
6. **Sample Diluent (PBS)**: 2 X 60 mL, dyed red, ready-to-use.
7. **Washing Solution (PBS + Tween 20)**: 60 mL as 10x concentrate. Dilute 1+9 with distilled water. If during the cold storage crystals precipitate, the concentrate should be warmed up at 37°C for 15 minutes.
8. Plastic bag to store unused microtiter strips.
9. Instruction Manual.

Instrumentation required but not provided

1. 50, 100 and 500 µL-micropipettes
2. ELISA reader (450 nm)
3. Centrifuge
4. Ultra-Turrax, mixer, vortex

Reagents required but not provided

1. Methanol
2. Double distilled water

Sample Preparation

Cereals / Nuts / Fruit

- Grind sample to pass through a 20 mesh sieve and thoroughly mix prior to sub-sampling.
- Suspend 20 g of sample in 100 mL of 70% methanol.
- Mix suspension for 5 minutes.
- Filter through Whatman #1 filter or alternatively centrifuge at a minimum of 3000g for 5 minutes.
- Dilute 500 µL of filtrate/supernatant with 500 µL of sample diluent and test the sample in the ELISA.

Dilution factor = 10

Beer / Gyle / Milk

- Dilute an adequate volume of sample diluent with 35% methanol.
- Defat milk if applicable. Therefore the milk has to be centrifuged for 15 min at 4°C and at least 2000g. Afterwards, the upper fat layer should be removed.
- Carbonized beer samples should be preliminarily degassed by moderate heating.
- Cloudy beers (such as beer brewed from wheat) / gyle should preliminarily be sterile-filtered.
- Dilute 100 µL of prepared liquid sample with 900 µL sample diluent/methanol dilution.

Dilution factor = 10

In case of samples having concentrations that are too high, an adequate volume of sample diluent is diluted with 35% methanol. The sample extracts have to be further diluted with this dilution.

ASSAY PROCEDURE

1. Prepare samples as described above.
2. Pipet 100 µL standards or prepared samples in duplicate into the appropriate wells of the microtiter plate.
3. Add 50 µL of zearalenone-peroxidase conjugate into each well.
4. Incubate for 30 minutes at room temperature.
5. Wash the plate three times as follows: Discard the contents of the wells (dump or aspirate). Pipet 300 µL of diluted washing solution into each well. After the third repetition empty the wells again and remove residual liquid by striking the plate against a paper towel. The wash procedure is critical. In-sufficient washing will result in poor precision and falsely elevated absorbencies.
6. Pipet 100 µL of substrate solution into each well.
7. Allow the reaction to develop in the dark (e.g. cup-board or drawer; the chromogen is light-sensitive) for 15 minutes at room temperature.
8. Stop enzyme reaction by adding 100 µL of stop solution (0.5 M H₂SO₄) into each well. The blue color will turn yellow upon addition.
9. After thorough mixing, measure absorbance at 450 nm (reference wavelength 620 nm), using an ELISA reader. The color is stable for 30 minutes.

RESULTS

1. Calculate the average optical density (OD 450 nm) for each set of reference standards or samples.
2. Construct a standard curve by plotting the mean optical density obtained for each reference standard against its concentration in ng/mL on semi-log graph paper with the optical density on the vertical (y) axis and the concentration on the horizontal (x) axis. Alternatively the evaluation can be carried out by software. In this case the 4-parameter method should be preferred.
3. The diluted samples must be further converted by the appropriate sample dilution factor for calculating the sample concentration in ppb. The factors for each sample matrix are listed in the sample preparation section.

Example:


A milk sample prepared as described above results in 0.6 ng/mL. The concentration of the sample is calculated as follows:

$$C_{\text{sample}} = 0.6 \text{ (ng/mL)} * 10 \text{ (ppb*ml/ng)} = 6 \text{ ppb}$$

Typical Standard Values

The following table contains an example for a typical standard curve. The binding is calculated as percent of the absorption of the 0 ng/mL standard. These values are only an example and should not be used instead of the standard curve which has to be measured in every new test.

Zearalenone (ng/mL)	(% binding of 0 ng/mL)
0	100
0.05	87
0.15	71
0.5	41
1.5	20
3.5	9

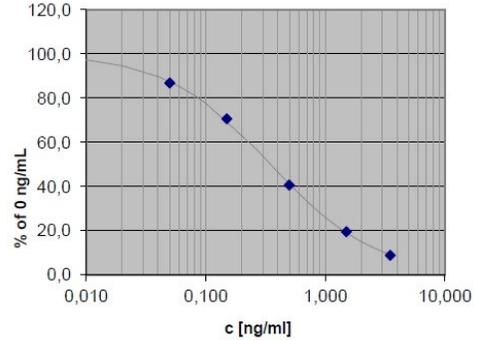

PERFORMANCE CHARACTERISTICS

Sensitivity

The limit of detection (LOD) of the Zearalenone test is 0.043 ng/mL.

The limit of quantification (LOQ) of the Zearalenone test is 0.130 ng/mL.

Validation experiments with common matrices result-ed in the following LODs and LOQs [ppb].

Matrix	LOD	LOQ
Wheat	0.7	1.2
Barley	0.9	2.0
Rye	0.9	1.9
Rice	0.7	1.7
Corn	0.5	1.4
Banana	1.5	3.5
Peanut	3.9	7.6
Milk	2.8	5.8
Beer	0.4	1.0

Recovery

Matrix	%
Wheat	86
Barley	103
Rye	105
Oats	110
Rice	93
Corn	101
Banana	86
Peanut	90
Milk	92
Beer	94

Linearity

The serial dilution of spiked samples (div. cereals, banana, peanut, milk and beer) resulted in a dilution linearity of 78-112%.

Precision

Intra-assay Precision	3-4%
Inter-assay Precision	3-13%

Cross-reactivity related to Zearalenone (=100%)

α-Zearalanol	4%
β-Zearalanol	26%
α-Zearalenol	6%
β-Zearalenol	107%

PRECAUTION

Full compliance of the following good laboratory practices (GLP) will determine the reliability of the results:

1. Prior to beginning the assay procedure, bring all reagents to room temperature (20-25°C).
2. All reagents should be mixed by gentle inversion or swirling prior to use. Do not induce foaming.
3. Once the assay has been started, all subsequent steps should be completed without interruption and within the recommended time limits.
4. Replace caps in all the reagents immediately after use. Do not interchange vial stoppers.
5. Use a separate disposable tip for each specimen to prevent cross-contamination.
6. All specimens and standards should be run at the same time, so that all conditions of testing are the same.
7. Do not mix components from different batches.
8. Do not use reagents after expiration date.
9. Check both precision and accuracy of the laboratory equipment used during the procedure (micropipets, ELISA reader etc.).

Health and safety instructions

1. Do not smoke or eat or drink or pipet by mouth in the laboratory.
2. Avoid contact of substrate and stop solution with skin and mucosa (possible irritation, burn or toxicity hazard). In case of contact, rinse the affected zone with plenty of water.
3. Handling and disposal of chemical products must be done according to good laboratory practices (GLP).

REFERENCES

1. Suzuki T, et al. (2007) – Sensitive detection of estrogenic mycotoxin zearalenone by open sandwich immunoassay. *Anal Sci*, 23(1):65-69
2. Sang-Ho C, et al. (2012) – Production of highly group-specific monoclonal antibody against zearalenone and its application in an enzyme-linked immunosorbent assay. *J Vet Sci*, 13(2):119-125
3. Wang YP, et al. (2006) – Development of ELISA-kit for quantitative analysis of zearalenone. *Wei Sheng Yan Jiu*, 35(2):221-224
4. Rashedi M, et al. (2012) – Zearalenone contamination in barley, corn, silage and wheat bran. *Toxic Ind Health*, 28(9):779-782
5. Gao Y, et al. (2012) – Preparation of highly specific anti-zearalenone antibodies by using the cationic protein conjugate and development of an indirect competitive enzyme-linked immunosorbent assay. *Analyst*, 137(1):229-236
6. Burmistrova NA, et al. (2009) – Application of a new anti-zearalenone monoclonal antibodies in different immunoassay formats. *Anal Bioanal Chem*, 395(5):1301-1307
7. Thongrussamee T, et al. (2008) – Monoclonal based enzyme-linked immunosorbent assay for the detection of zearalenone in cereals. *Food Add Contam Part*, 25(8):997-1006
8. Bennett GA, et al. (1994) – Enzyme-linked immunosorbent assay for detection of zearalenone in corn, wheat and pig feed: collaborative study. *J AOAC Int*, 77(6):1500-1509

9. Vulich A, et al. (2012) – Analysis of naturally occurring zearalenone in feeding stuffs and urine of farm animals in Croatia. *J Immunoassay Im-munochem*, 33(4):369-376
10. Hervás M, et al. (2009) – Electrochemical micro-fluidic chips coupled to magnetic bead-based ELISA to control allowable levels of zearalenone in baby foods using simplified calibration. *Ana-lyst*, 134(12):2405-2411

<p>ISO 13485 ISO 9001</p> 	
 <p>Diagnostic Automation/Cortez Diagnostics, Inc. 21250 Califa Street, Suite 102 and 116, Woodland Hills, California 91367 USA</p>	
Date Adopted	2017-9-20
REF 5148-8	AccuDiag™- Zearalenone ELISA
Revision Date: 2017-3-31	